*Well
Store*®

Estetico & Benessere

CE

HSD

ITA - ENG - FRA - SPA - DEU

ITALIANO	Programmi per l'uso estetico & benessere	PAGE 3-5
ENGLISH	Programs for Aesthetic & Wellness use	PAGE 8-10
FRANÇAIS	Programmes à usage esthétique & bien-être	PAGE 11-14
ESPAÑOL	Programas de uso estético & bienestar	PAGE 15-18
DEUTSCH	Programme für Beauty & Wellness	PAGE 19-22

----- ITALIANO -----

Programmi Estetici-Benessere

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PROGRAMMA 98

----- ENGLISH -----

Programs for Aesthetic-Wellness use

Information

In programs where the area of the body to be treated allows it, the use of the ultrasonic handpiece can be replaced with double-head diffusers.

PROGRAMS

Application notes to use the device correctly, knowing which program to choose and the action it performs, thus fighting your imperfection as effectively as possible.

Bump

ACTION The ultrasound irradiation in this action causes a massage of considerable intensity, which by means of the oscillations penetrates deeper into the epidermis.

WHEN It should be used in all those advanced stages of localized cellulite, especially in the case of soft or flaccid cellulite.

Modeling

ACTION This therapy is widely used in cases of localized lipodi-strophy by breaking the fibronic tissue present which hinders both blood and lymphatic circulation.

WHEN Benefiting from its characteristic in shaping and stimulating it is particularly indicated in the case of hard or solid cellulite.

Maintenance

ACTION Its stimulating effect on local blood circulation improves skin trophism and cellular oxygenation, releasing excess fluids and toxins from the subcutaneous connective tissue.

WHEN It is indicated as maintaining the cycle of sessions carried out, in order not to lose the results achieved.

Facial scrub

ACTION It reactivates the circulation and promotes cell turnover: the skin is purified in depth and prepared to better absorb the active ingredients of moisturizing and anti-wrinkle creams. The exfoliating action also has beneficial effects on blackheads which are reduced or even disappear with ultrasonic waves. The treatment should be performed once a day.

Face massage

ACTION It increases the elasticity and tone of the tissues by reducing the possible formation of post-acne scars and increases the effect of anti-acne products. The treatment should be performed once a day.

Smoothing for the face

ACTION It smoothes and relaxes (expression) wrinkles and if used preventively it can delay skin aging or prevent it from worsening. Thanks to its thermal effects, ultrasound heats the tissues and promotes biochemical and metabolic processes. The treatment should be performed once a day.

Smoothing

ACTION It recalls the clean plasma throughout the thickness of the dermis and restores the right physiological dimension to the skin connective tissue, improving both the elasticity and the superficial irregularity of the epidermis.

WHEN It is indicated in all those cellulite blemishes, in which the skin surface has irregularities, such as the classic "Orange peel skin" and "Rough Skin".

Draining

ACTION Its beneficial action facilitates the drainage of the lymphatic connectors and of the respective attached glandular stations, facilitating the reabsorption of liquids.

WHEN Benefiting from its draining action, it is recommended in the case of edematous cellulitis, as liquids are the main cause, also recommended against orange peel skin.

Reactivating

ACTION Stimulates and reactivates blood and lymphatic circulation for skin improvement, reactivating cell metabolism, favoring the exchange of liquids between cells.

WHEN It is especially indicated in the initial stages of cellulite, as it stimulates and reactivates the dormant metabolism.

Carrier

ACTION It favors the penetration of creams, gels or other products into the biological tissue, favoring their absorption and therefore their effect, also indicated for depilatory treatments.

PROGRAM

OUTER ARMS

ACTION

SHOCK	01
MODELING	02

PROGRAM

INNER ARMS

ACTION

SHOCK	03
MODELING	04

PROGRAM

HIPS AND BACK

ACTION

SHOCK	05
MODELING	06

PROGRAM


ABDOMIN. AREA

ACTION


SHOCK	07
MODELING	08

PROGRAM	
CULOTT. CHEVAL	
ACTION	
SHOCK	09
MODELING	10

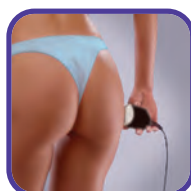PROGRAM	
GLUTEUS MUSCLE	
ACTION	
SHOCK	11
MODELING	12

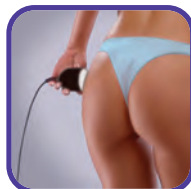PROGRAM	
OUTER THIGHS	
ACTION	
SHOCK	13
MODELING	14

PROGRAM	
INNER THIGHS	
ACTION	
SHOCK	15
MODELING	16

PROGRAM	
KNEES	
ACTION	
SHOCK	17
MODELING	18

PROGRAM	
INNER KNEE	
ACTION	
SHOCK	19
MODELING	20

PROGRAM	
CALVES	
ACTION	
SHOCK	21
MODELING	22


PROGRAM	
CELLULITE	
ACTION	
DRAINING	23
SMOOTHING	24

PROGRAM	
ROUGH SKIN	
ACTION	
REACTIVATING	25
SMOOTHING	26

PROGRAM	
MAN	
ACTION	
SHOCK	27
MODELING	28
MAINTENANCE	29

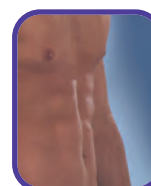PROGRAM	
VEHICULATING	
ACTION	
HIGH	30
MEDIUM	31
LOW	32

PROGRAM	
FACE	
ACTION	
EXFOLIATING	33
MASSAGE	34
LEVIGATING	35

PROGRAM	
SCALP	
ACTION	
STIMULATING	36
CLEANING	37
ACTIVATING	38

PROGRAM	
IMMERSION	
ACTION	
1 MHz	39
3 MHz	40

PROGRAM	
ANTI-AGE ACTION	
ACTION	
ANTI-AGE	41

PROGRAM	
VARICOSE VEINS TREATMENT	
ACTION	
VARICOSE VEINS	47

PROGRAM	
SKIN SPOTS	
ACTION	
SKIN SPOTS	42

PROGRAM	
ACNE AND SPOTS	
ACTION	
ACNE AND SPOTS	48

PROGRAM	
SKIN TISSUE REGENERATION	
ACTION	
SKIN TISSUE REGENERATION	43


PROGRAM	
CELLULITE ABDOMEN	
ACTION	
CELLULITE ABDOMEN	49

PROGRAM	
WRINKLE TREATMENT	
ACTION	
WRINKLE TREATMENT	44

PROGRAM	
CELLULITE LOWER LIMBS	
ACTION	
CELLULITE LOWER LIMBS	50

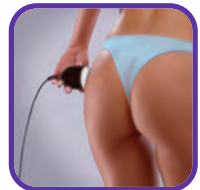PROGRAM	
COUPEROSE TREATMENT	
ACTION	
COUPEROSE TREATMENT	45


PROGRAM	
CELLULITE UPPER LIMBS	
ACTION	
CELLULITE UPPER LIMBS	51

PROGRAM	
STRETCH MARKS TREATMENT	
ACTION	
STRE. MARKS TREATMENT	46

PROGRAM	
CELLULITE BUTTOCKS	
ACTION	
CELLULITE BUTTOCKS	52

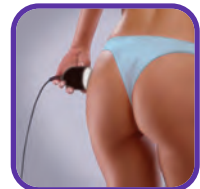PROBE TESTING

12 minute program to be used only and exclusively to check the correct functioning of the treatment handpiece

PROGRAM	98
---------	----

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